

Recommended Steps for Putting On and Taking Off Personal Protective Equipment

This resource is an excerpt from the [Routine Practices and Additional Precautions, In All Health Care Settings \(Appendix L\)](#) and was reformatted for ease of use.

For more information please contact Public Health Ontario's Infection Prevention and Control Department at ipac@oahpp.ca or visit www.publichealthontario.ca

[Images developed by Kevin Rostant. Some images adapted from Northwestern Ontario Infection Control Network – NWOICN]

This is an excerpt from [Routine Practices and Additional Precautions In All Health Care Settings \(Appendix L\)](#)


PUTTING ON PPE

1. Perform Hand Hygiene




2. Put on Gown

- Tie neck and waist ties securely




5. Put on Gloves

- Put on gloves, taking care not to tear or puncture glove
- If a gown is worn, the glove fits over the gown's cuff



3. Put on Mask/N95 Respirator

- Place mask over nose and under chin
- Secure ties, loops or straps
- Mould metal piece to your nose bridge
- For respirators, perform a seal-check



4. Put on Protective Eyewear

- Put on eye protection and adjust to fit
- Face shield should fit over brow




Routine Practices and Additional Precautions (November 2012)

This is an excerpt from [Routine Practices and Additional Precautions In All Health Care Settings \(Appendix L\)](#)

TAKING OFF PPE

1. Remove Gloves

- Remove gloves using a glove-to-glove/skin-to-skin technique
- Grasp outside edge near the wrist and peel away, rolling the glove inside-out
- Reach under the second glove and peel away
- Discard immediately into waste receptacle

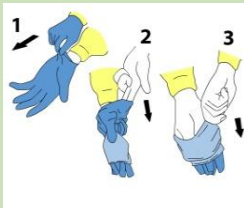

2. Remove Gown

- Remove gown in a manner that prevents contamination of clothing or skin
- Starting at the neck ties, the outer, 'contaminated', side of the gown is pulled forward and turned inward, rolled off the arms into a bundle, then discarded immediately in a manner that minimizes air disturbance



6. Perform Hand Hygiene




3. Perform Hand Hygiene



5. Remove Mask/N95 Respirator

- Ties/ear loops/straps are considered 'clean' and may be touched with hands
- The front of the mask/respirator is considered to be contaminated
- Untie bottom tie then top tie, or grasp straps or ear loops
- Pull forward off the head, bending forward to allow mask/respirator to fall away from the face
- Discard immediately into waste receptacle

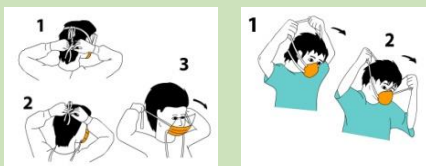

4. Remove Eye Protection

- Arms of goggles and headband of face shields are considered to be 'clean' and may be touched with the hands
- The front of goggles/face shield is considered to be contaminated
- Remove eye protection by handling ear loops, sides or back only
- Discard into waste receptacle or into appropriate container to be sent for reprocessing
- Personally-owned eyewear may be cleaned by the individual after each use

